

Specifications

Power Supply	100~240VAC, 50/60Hz, 200VA
Flow Rate	0.1~10ml/s
Pressure Limit	325PSI
Auto Filling Rate	8ml/s
Air Purge Rate	A: 1.0~9.9ml/s B: 2.5~9.9ml/s
Piston Rod Speed	9.9ml/s
Injection Delay	1~300s
Scan Delay	1~300s
Pause Phase	1~900s
Phases	6
Injection Programs	100
Injection Records	100
Volume	A: 1~65ml
	B: 1~115ml
Multiple Language User Interface	



DC Power

Medica

High Quality Disposables



+786-9201598

sales@globalmedica.us

- 65ml Syringe for Contrast
- · 115ml Syringe for Saline





MR Injector

Leading Contrast Media Injectors Forward



MR Injector


Introducing an all-new ImaStar MR Injector, offering an extended list of product features, reliable performance, and a safe patient experience.



- Ease of Operation and Clinical Efficiency
- · Fill contrast and saline with a single press of a button
- Optional IV pole

Independent Dual Motor Drive Systems

- Dual motor drive capability for time saving loading and filling of syringes
- · Allows independent injection of contrast and saline

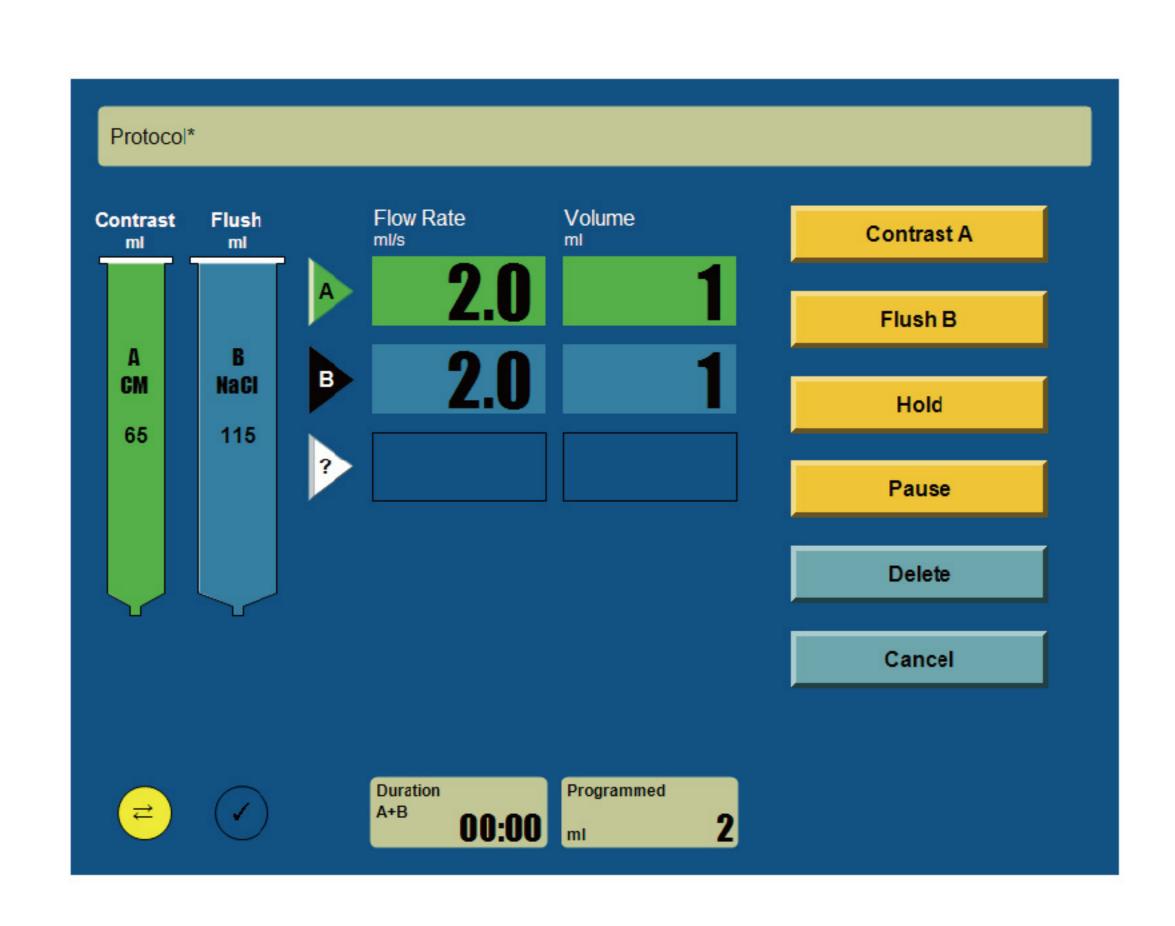



Color Touch-Screen

- · 12.1"color touch-screens for easy operation
- Clear and concise information display to ensure user operational confidence


Injection Record Storage

· Stores up to 100 injection history records



Real-time Pressure Curve


 Display provides real-time pressure graph for injection monitoring

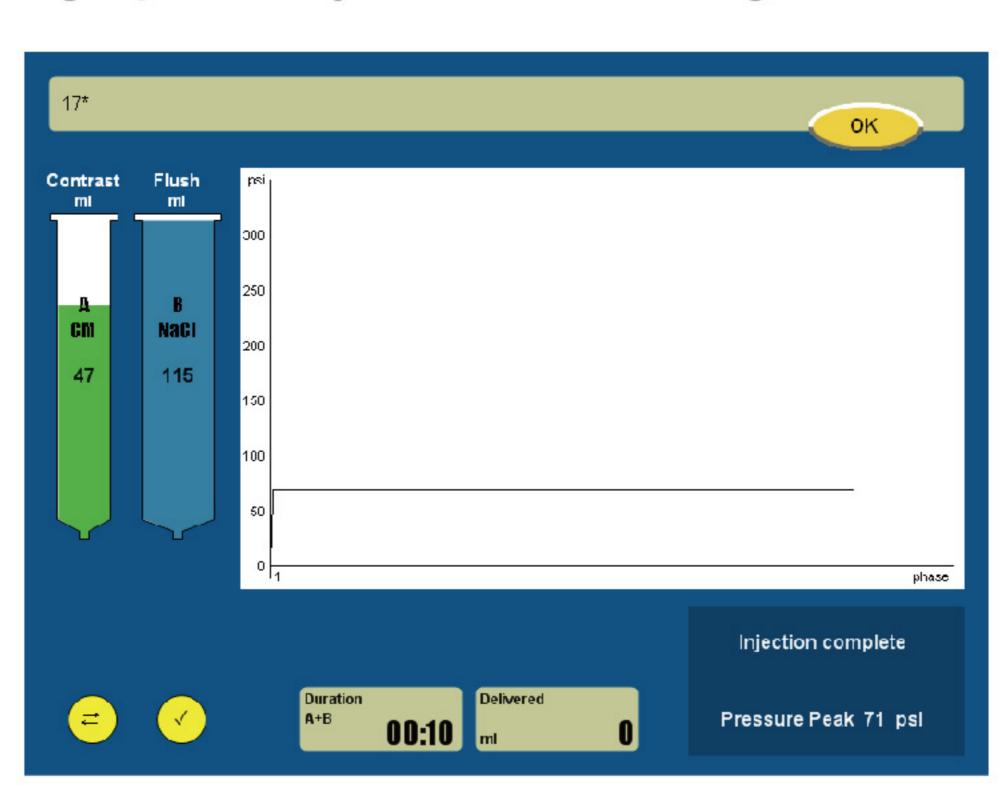

Configuration Options


A Variety of

- Stores up to 100 customized injection programs
- Program storage capacity can be increased based on user needs

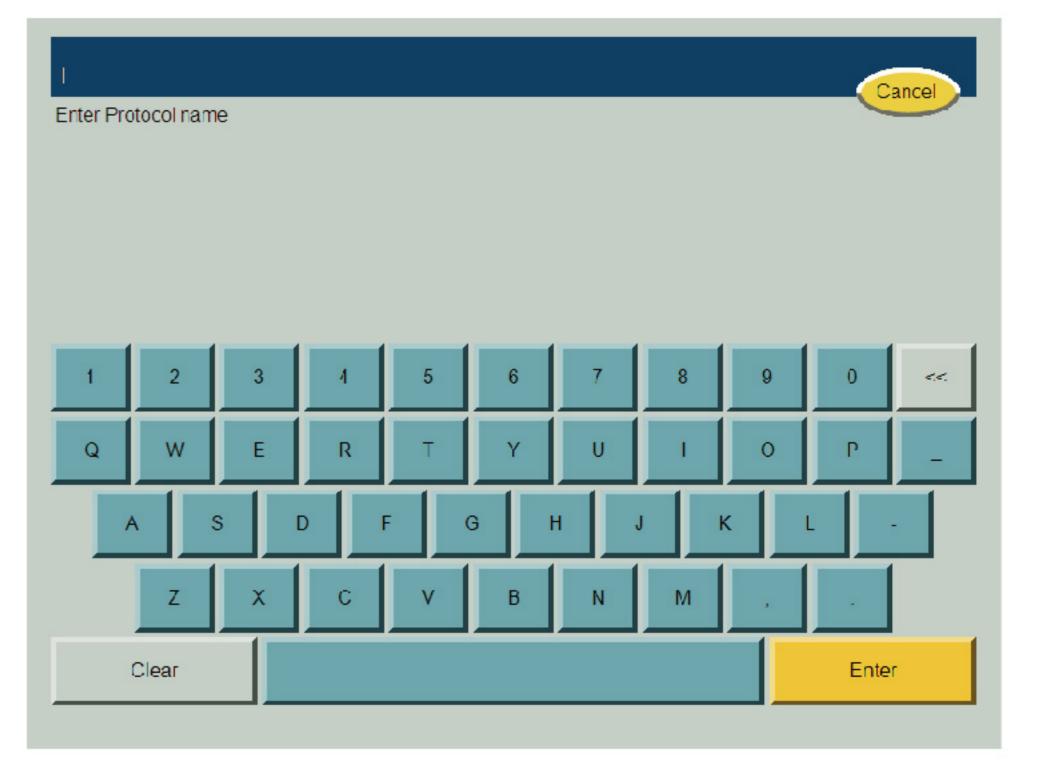

KVO Function

- Independent KVO (Keep Vein Open) ensures efficient examination procedure
- Flexible KVO interval time settings
- 115ml saline syringe capacity allows for longer KVO

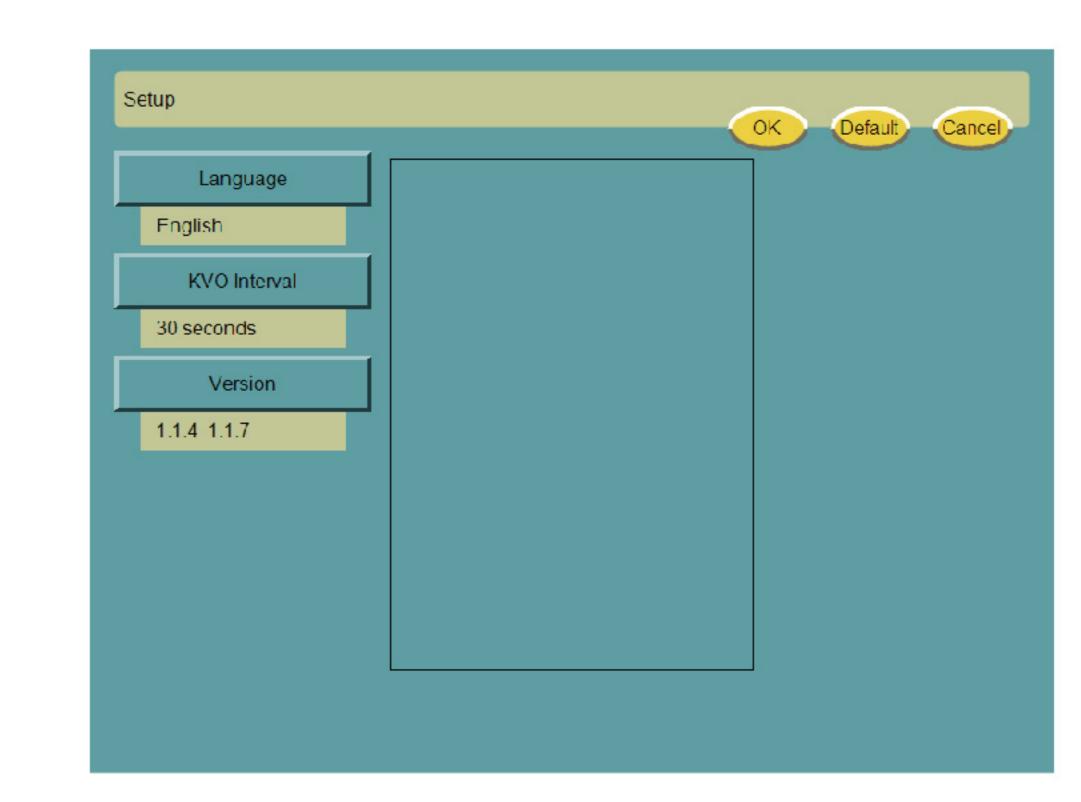

Safety

- To ensure the safety of the injection, there is a double confirmation feature with separate prompts to confirm air purging and injection ready status before any injection
- · Automatically slows down the injection when approaching to pressure limit
- · Injection stops when pressure limit is reached. Alerts with both light and sound
- Waterproof design inside the head to minimize damage from contrast media and saline leakage